HYUNDAI

OVET 111

NÁVOD K POUŽITÍ / NÁVOD NA POUŽITIE INSTRUKCJA OBSŁUGI / INSTRUCTION MANUAL HASZNÁLATÍ ÚTMUTATÓ

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Producent: Hyundai Corporation, Seoul, Korea.

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Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.


GENERAL SAFETY INSTRUCTIONS

- The appliance is designed exclusively for private use and for the envisaged purpose.
- It is not intended for use in health-care facilities and for commercial use!
- Keep it away from sources of heat, direct sunlight, humidity (never dip it into any liquid) and sharp edges. Do not use the appliance with wet hands.
- If you should leave the workplace, always switch the machine off.
- Do not try to repair the appliance on your own. Always contact an authorized technician.
- Pay careful attention to the following “Special Safety Instructions for this unit“.

SPECIAL SAFETY INSTRUCTIONS FOR THIS UNIT

- Operate the appliance only in accordance with the specifications on the rating plate.
- Do not put the appliance on the wet places.
- Clean the scale with a slightly damp cloth.
Do not immerse the scale in the water or use chemical / abrasive agents.
- Always use the scale on a hard, flat surface. **DO NOT** use on carpet
- **DO NOT** use the scale on wet and slippery floor. **DO NOT** step on the edge of the scale.
The scale platform is slippery when wet. Keep it dry!
- Always weigh yourself with the same scales places on the same surface.
Uneven floors may affect the reading.
- Weigh yourself without clothing and footwear, before meals and always at the same time of day.
- Do not use this scale if you are using one of the following devices: pace maker, electrocardiograph and/or other medical device implanted in the body or device used as support.
- The body fat & hydration are for reference only. Please always consult your doctor for advice of medication or diet in order to achieve ideal body fat level.
- To make sure accuracy of testing, please **DO NOT** use the scale in environment with strong electro-magnetic field.
- **DO NOT** use the scale on wet and slippery floor. **DO NOT** step on the edge of the scale.
The scale platform is slippery when wet. Keep it dry!
- **DO NOT** use this scale on a place with intense shake. Keep stand still on scale when testing.
- Do not drop the device, or jump on the device when using. Avoid strong impact and sudden shake as it is a precision measuring device. The construction is delicate and sensitive to all unwanted external impact.
- Keep in mind that the data read from this scale are approximate and inaccurate in comparison to results of actual medicinal analyses. The exact determination of biometrical characteristics can be done only by a doctor via medical methods.
- Keep batteries and scale out of reach of children and incapacitated people. A person, who would swallow the batteries, must seek medical help immediately.
- Do not use during pregnancy.
- The appliance can be used by people (including children) whose physical, sensual or mental inability or lack of experience and knowledge prevents safe use of the appliance if they are supervised.
- Children should be supervised to ensure that they do not play with the appliance.

DESCRIPTION OF THE CONTROLS

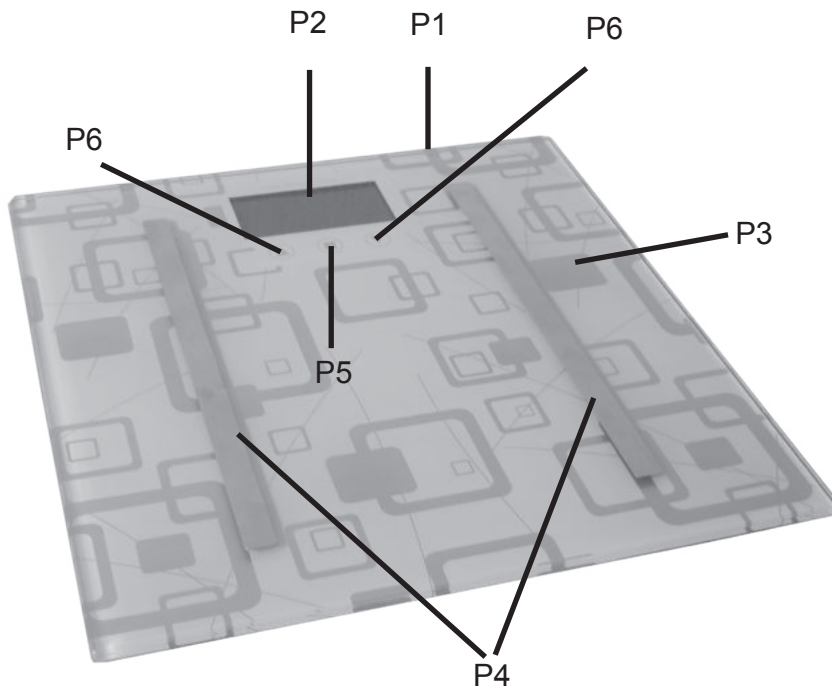- P1. Unit of scale
- P2. Display
- P3. Glass surface
- P4. Electrodes
- P5. SET / ON button
- P6. UP and DOWN buttons for setting

LCD SEGMENTUM / ICONS ILLUSTRATIONS:

P1 ...P0 number of memory

man

woman

musculature contents (%)

bone weight contents (%)

fat body fat contents (%)

Underfat

Healthy

Overfat

Obese

TBW body water contents (%)

KCAL optimum energy consumption

SPECIFICATIONS:

Max. Weighing Capacity: 180 kg, 28 st (stone), 396 lb (libra)

Weight Measuring Division: 0,1 kg, 0,2 lb, 1/4 st

Water division: 0,1 %

Body fat division: 0,1 %

Age range for body fat: 6 - 100 years

Height range: 80 - 220 cm / 3ft 3,5in - 7 ft 2,5 in

Height unit: cm / ft-in

INSTRUCTIONS FOR USE

BEFORE USE

- Remove the insulation tab from the battery contacts before use.
- Activate the scale and wait until it shows "0.0". Choose a measuring unit (KG/ST/LB) with UP and DOWN buttons.
- Place the scales on a firm, flat and even surface to ensure accuracy.

WEIGHT MEASURING

You can use this scale only for weight measuring.

This is basic mode without body fat measuring.

- 1) Gently step on the scale, then the scale will be automatically turned on. Stand evenly on the scale without moving and wait until your weight shown on display is stable and locked. When you step off the scales it will automatically switch off after 8 seconds.
- 2) If no weighing is carried out while 0.0 is displayed the scales will automatically switch off after 8 seconds.

BODY FAT MEASUREMENT

- 1) About Body Fat Percentage

Body fat percentage indicates body fat proportion in body, which is a criteria to judge one is overweight or not.

- 2) Testing Principle

Impedance of fat is higher than any other body tissues such as muscle, when electrical signal of a certain frequency goes through body (BIA principle). Based on this principle, body fat percentage can be tested by sending an electrical signal through body.

BMI

The technology used in these scales, Bio-electrical Impedance Analysis (BIA) is more technologically advanced than the order Body Mass Index (BMI) method where fat % is calculated using only height and weight. BIA technology uses 5 parameters to get a more accurate resolution for the analysis (weight and biological resistance which are automatically measured by the scales together with the three personal parameters entered manually - sex, height and age). This index sets mutual ratio between the body weight and height. It can be used to determine whether the person suffers from underweight, overweight or obesity or whether his/her health condition is normal. Also with regards to heart and vascular diseases, body mass index is a relatively good risk index, particularly in healthy people. However, the index is for reference only. Experts recommend supplementing the calculation by at least measuring the waist, which makes determining disposition of the body fat more precise. BMI will not provide very reliable information for small children, elderly people and active sportsmen who have a lot of muscle mass as their results can be biased.

Optimum energy consumption CAL

It represents optimum amount of energy consumption necessary for maintaining the basic vital functions and balanced metabolism. The data on the display show the optimum daily consumption of energy that you give out at your weight, height, age and gender. You should also take in this amount of energy to keep your metabolism balanced. If you intend to reduce your weight, your daily energy intake should be a bit lower. However, remember to keep carbohydrates, fats and proteins balanced in your diet.

239 cal = 1 kJ

Important:

During body fat testing, please remove shoes and socks, etc., and clear your feet, in order to get the most accurate result.

To ensure accuracy, please test at the same time of the day

Test result may be inaccurate when the users are:

- after intense sport activity, overfeeding, or excessive dehydration,
- children under 10, or adults over 80,
- body building athletes, or other occupational athletes,
- pregnant women; Dialysis or edema patient,
- do not use the fat analysis feature on this scale if you are using one of the following devices: pace maker, electrocardiograph and/or other medical device implanted in the body or device used as support,
- excessive drinking, eating, exercising, health conditions, used medicaments, women's menstruation cycle etc. influence the results of the analysis,
- figyeljen arra, hogy mérés közben ne legyen kontaktus a talpak, vádlik, vagy combok között. Ellenkező esetben lehetséges, hogy a mérés nem lesz pontos.

In order to ensure accurate results users should enter the required personal parameters; age, sex and height. These scale have the facility to store the personal parameters of 10 people.

SETTING AND SAVING PARAMETERS

Press to turn ON.

Press or to switch between Body Parameter P0 - P9. Press to confirm.

Press or to select male or female. Press to confirm.

Press or to indicate height. Press to confirm.

Press or to indicate age. Press to confirm.

The values are flashed and recorded.

MEASURING

- 1) Tap on the scales platform to switch on the product
- 2) Press P1 will flash. Press or to switch between Body Parameter P0 - P9.
- 3) Step onto the platform when the LCD display reads 0.0.
Important: For correct measurements and scale functions step on the scale barefoot.
- 4) First the weight is recorded, then „ o o o o “ is displayed, indicating that the scales are analyzing. Please do not alight.
- 5) Once finished, the weight, body fat, body muscle and bone weight will be displayed alternately for 2 times before the scales automatically turn off.
- 6) After this the scale is going turn off.
- 7) If no weighing is carried out while 0.0 is displayed the scales will automatically switch off after 8 seconds.

Compare the body fat reading against the <Body Fat Ratio Chart> to confirm the body status.

Body Fat Ratio Chart

| Age | Female (%) | | | | Male (%) | | | |
|----------|------------|---------|---------|-------|----------|---------|---------|-------|
| | Low | Healthy | High | Obese | Low | Healthy | High | Obese |
| 6 - 39 | < 21 | 23 - 33 | 33 - 39 | > 39 | < 8 | 8 - 19 | 19 - 25 | > 25 |
| 40 - 30 | < 23 | 23 - 35 | 35 - 40 | > 40 | < 11 | 11 - 22 | 22 - 28 | > 28 |
| 60 - 100 | < 24 | 24 - 36 | 36 - 42 | > 42 | < 13 | 13 - 25 | 25 - 30 | > 30 |

Water Ratio Chart & Muscle Ratio Chart

| Age | Female % water | | | Body Muscle | Male % water | | | Body Muscle |
|----------|----------------|---------|------|-------------|--------------|---------|------|-------------|
| | Low | Healthy | High | | Low | Healthy | high | |
| 6 - 15 | < 57 | 57 - 67 | > 67 | > 34 | < 58 | 58 - 72 | > 72 | > 40 |
| 16 - 30 | < 47 | 47 - 57 | > 57 | > 34 | < 53 | 53 - 67 | > 67 | > 40 |
| 31 - 60 | < 42 | 42 - 52 | > 52 | > 34 | < 57 | 47 - 61 | > 61 | > 40 |
| 61 - 100 | < 37 | 37 - 47 | > 47 | > 34 | < 42 | 42 - 56 | > 56 | > 40 |

Bone Weight Chart

| Weight | Female | | | Male | | |
|------------------|----------|----------|----------|----------|----------|----------|
| | < 45 kg | 45-60 kg | > 60 kg | < 60 kg | 60-75 kg | > 75 kg |
| Ref. bone weight | ~ 3.0 kg | ~ 4.2 kg | ~ 6.5 kg | ~ 4.5 kg | ~ 6.0 kg | ~ 7.5 kg |

KCAL Chart

| Věk | Female | | Male | |
|----------|-------------|------|-------------|------|
| | Weight (kg) | KCAL | Weight (kg) | KCAL |
| 6 - 17 | 50 | 1265 | 60 | 1620 |
| 18 - 29 | 55 | 1298 | 65 | 1560 |
| 30 - 49 | 60 | 1416 | 70 | 1680 |
| 50 - 69 | 60 | 1302 | 70 | 1561 |
| 70 - 100 | 50 | 1035 | 65 | 1397 |

BMI Chart

| BMI | Underfat | Healthy | Overfat | Obese |
|-----|----------|-----------|---------|-------|
| | < 18.5 | 18.5 - 27 | 30 - 39 | > 40 |

Due to the natural differences between people this chart should be used for reference purposes only.

ERROR AND OTHER LCD MESSAGE

LCD may show occasionally error message 'Err',. Simply turn off the unit for a brief moment, or remove and re-install battery as a total re-set. Then turn on again to resume operation.

- 1) When the LCD screen displays „LO“, you should install a new battery.
- 2) When the LCD screen displays „Err“, it indicates that the maximum weight is exceeded.
- 3) When the LCD screen displays „Err 2“, it shows an error in the body fat percentage. Fat percentage exceeds the maximum capacity of 50 %, fat percentage is under 3 %.

If signs appear, such as if the display is flashing or if there is a faulty sign, try to move the scale away from the influence of a source of interference or ensure switching off the source for the time of using the scale.

CLEANING AND MAINTENANCE

- 1) Use a dry, soft fabric to clean the device. Or a lightly damp soft fabric with mild detergent to wipe clean the device. Do not use strong chemical solution such as alcohol, benzene for cleaning.
- 2) It could be slippery when the unit platform surface is wet, or semi-damp. Ensure the weighing platform surface is dry before using for safety reason.

TROUBLESHOOTING

The device is not working

- Check the right position of battery and the polarity
- When warning indicator shows „Lo“ replace the battery.

TECHNICAL SPECIFICATIONS

- Bathroom body analysis scale
- Glass weighing surface
- Maximum weight capacity up to 180 kg
- Measurement accuracy 100 g
- Measurement units kg/lb/st
- 6 mm safety glass
- Weighs immediately - just step on the scales - it will automatically turn on, ready to weigh
- Automatic shut-off, automatic reset
- Overload indicator, low battery indicator
- „Using 5 parameters for measuring – height, weight, age, sex and biological resistant“
- It measures body fat, body water, muscles, bones, KCAL and BMI
- Data of 10 persons (gender, age, height) can be saved for analysis
- Sensors for higher weighing accuracy
- Large display
- Touch control
- Operation on battery type CR2032 (2 pc)
- (Included)
- Protection class III, CE Certification
- Color: white + silver
- Size: 30 x 2.2 x 30 cm
- Weight (netto): 1.7 kg

LEGISLATIVE & ECOLOGY

The product has received ES declaration of conformity according to Act as amended.

The product matches the requirements of the below statutory order as amended.

– Committee Regulation No. 2004/108/ES as amended, setting the technical requirements of products in terms of their electromagnetic compatibility.

The manufacturer reserves the right to make insignificant changes to the standard design that do not affect the function.

WARNING: Do not use this product near water, in wet areas to avoid fire or injury of electric current. Always turn off the product when you don't use it or before a revision. There aren't any parts in this appliance which are reparable by consumer. Always appeal to a qualified authorized service. The product is under a dangerous tention.

Old electrical appliances, used batteries and accumulators disposal


This symbol appearing on the product, on the product accessories or on the product packing means that the product must not be disposed as household waste. When the product/ battery durability is over, please, deliver the product or battery (if it is enclosed) to the respective collection point, where the electrical appliances or batteries will be recycled. The places, where the used electrical appliances are collected, exist in the European Union and in other European countries as well. By proper disposal of the product you can prevent possible negative impact on environment and human health, which might otherwise occur as a consequence of improper manipulation with the product or battery/ accumulator. Recycling of materials contributes to protection of natural resources. Therefore, please, do not throw the old electrical appliances and batteries/ accumulators in the household waste. Information, where it is possible to leave the old electrical appliances for free, is provided at your local authority, at the store where you have bought the product.


Information, where you can leave the batteries and accumulators for free, is provided to you at the store, at your local authority.

HOUSEHOLD USE ONLY DO NOT IMMERSE IN WATER OR OTHER LIQUIDS TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES, OR PLAYPENS. THIS BAG IS NOT A TOY.

 NOTICE

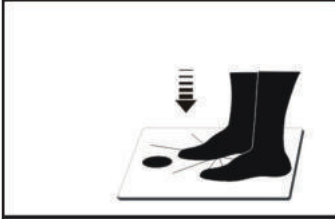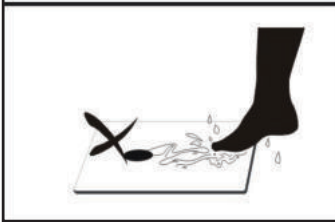AUTO STEP ON – just step on the scale and stand still, the scale will automatically turn ON and start taking measurements.

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Licensed by Hyundai Corporation, Korea.