

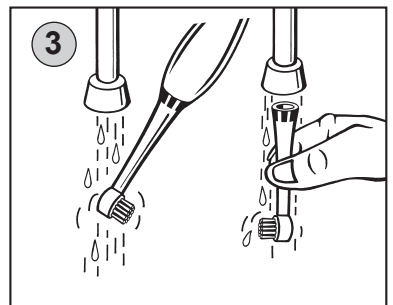
<i>Bateriový zubní kartáček • NÁVOD K OBSLUZE</i>	CZ 3-5
<i>Batériová zubná kefka • NÁVOD NA OBSLUHU</i>	SK 6-8
<i>Battery toothbrush • INSTRUCTIONS FOR USE</i>	GB 9-11
<i>Elektromos fogkefe • HASZNÁLATI UTASÍTÁS</i>	H 12-14
<i>Szczoteczka do zębów na baterie • INSTRUKCJA OBSŁUGI</i>	PL 15-17
<i>Зубная щётка • ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ</i>	RU 18-21

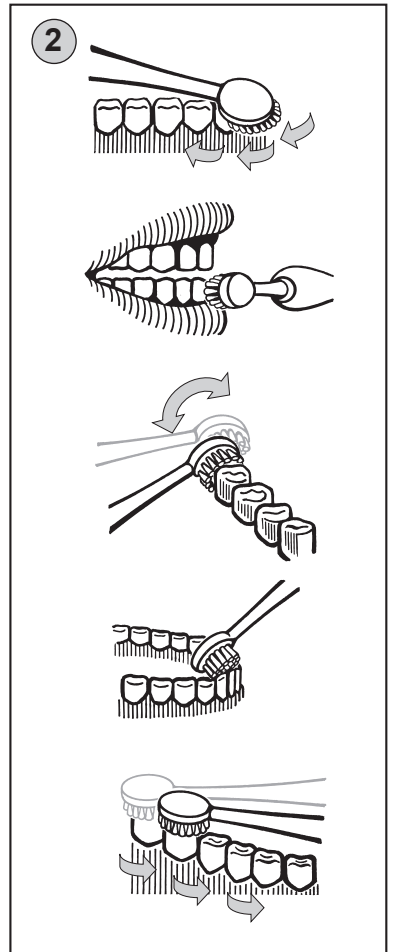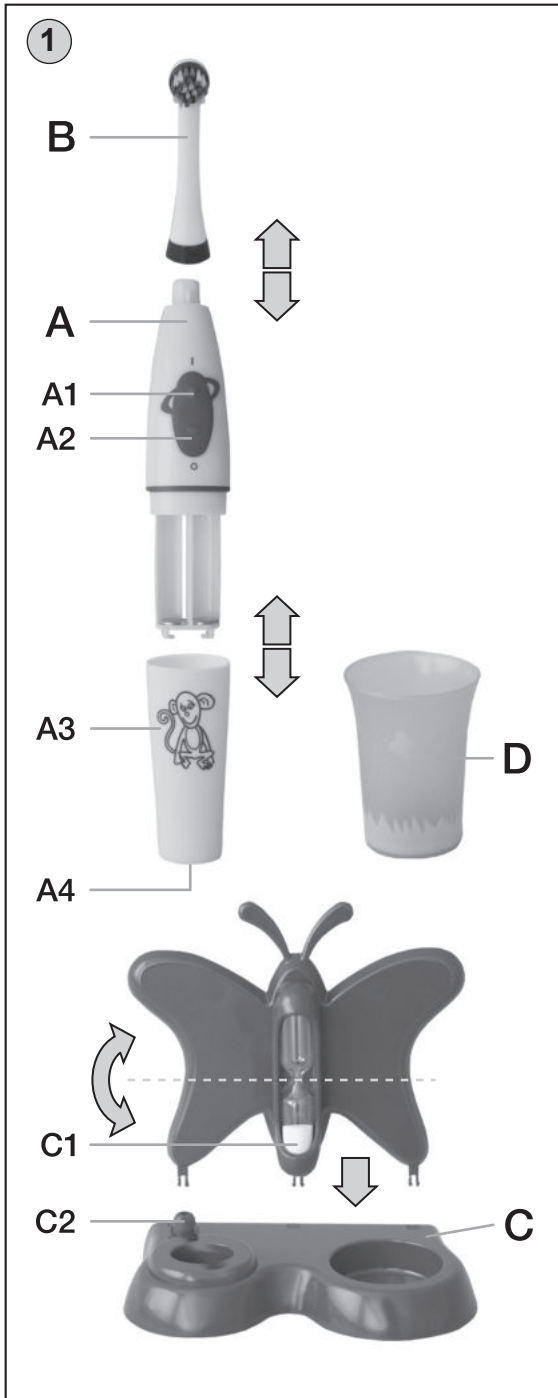eta 1294

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eta 1294

USER MANUAL

This device is a mechanical remover of tooth plaque which is often used and which will help your children improve their oral hygiene. Tooth brush represents the latest technology contained in a effective and ergonomic device which is an ideal tool for teeth care.

I. SAFETY PRECAUTIONS

- Before first use, please read this manual, look at the figures and keep it for further reference.
- Never use the device if it's not working properly, has fallen on the ground or has been damaged.
- The toothbrushes may be used by children and persons with decreased physical or mental abilities or lack of experience and knowledge, if they are under supervision or have been instructed on using the appliance in a safe way and understand potential dangers. Cleaning and maintenance carried out by the user may not be done by children without supervision. Children must not play with the appliance.
- Never immerse the propulsion unit in water and other liquids!
- Use the appliance only with accessory (replacement brush heads) intended for this type of device.
- Do not put the brush away on hot heat sources (e.g. **stove, radiator, cooker**)
- If the device has been stored at lower temperatures, acclimatize it first which will maintain the mechanical compactness of the plastic parts.
- Do not put the brush into any bodily cavities (except the oral cavity).
- Clean the brush after use (see par. **IV. MAINTANENCE**) and store it at a dry place.
- Use only the same type of batteries and do not combine used batteries with new ones.
- If you are not going to use the appliance for a longer period of time, take out the batteries.
- If the batteries are leaking, take them out immediately, otherwise they may damage the toothbrush
- Dispose of the old batteries appropriately (see paragraph **V. ECOLOGY**).
- Do not use the appliance for any other use than the one it is intended for and described!
- The manufacturer is not liable for any damage caused by incorrect use of the device (e.g. injury of oral cavity, gums, teeth, damaging the device) and is not liable for the warranties in case of breaching the precautions mentioned above.

II. Features of the device and accessories (pic. 1)

A - propulsion unit

A1 – button I (switch)

A2 – button 0 (switch)

A3 - battery cover

A4 – battery cover arrest control

B - brush head

C - stand

C1 – timer

C2 – brush head holder

D - cup

III. INSTRUCTIONS FOR USE

Remove all packaging material and remove the toothbrush with accessories. Rinse the parts that will come in touch with the oral cavity in clean hot water before first use.

Notes

For hygiene reasons we recommend every user using their own brush head.

You may use any quality tooth paste for cleaning your teeth. We do not, however, recommend using special tooth pastes with whitening effects. In order to prevent spraying, first put the brush on the teeth and then turn on the propulsion unit.

If you have undergone any surgical intervention in the mouth cavity, consult your dentist before using the toothbrush.

During the first few days of using the toothbrush, your gums may bleed slightly which, however, usually stops by itself after a few days. If the bleeding lasts longer than two weeks, please consult your dentist.

Inserting batteries

Release the battery cover **A3** on the bottom of the propulsion unit **A** by turning the control **A4** clock-wise into a position ● and remove the cover. Insert 2 pcs. of AAA batteries, observing the +/- signs (see symbol). Replace the cover and arrest it by turning the control into the position ☰.

Use

Slide the brush head **B** slightly on the propulsion unit **A** so that its protrusions fit onto the protrusions of the propulsion unit. Rinse your mouth thoroughly before cleaning your teeth, put toothpaste on the brush head and moisten it. Position the brush against the teeth so that it covers the surface of the teeth as well as the gum line. Turn the brush on by pressing the **I** button.

Guide the brush slowly along the teeth and grinders, follow the curve of the gums and the shape of the individual teeth. Hold the brush always on one spot for several seconds before you move over to a next tooth. Make sure you have cleaned the whole surface of the teeth (the inner and outer part of the occlusion area), as well as the line around the gums (pic. 2). Do not press on the brush, leave it to do all the work by itself. Acquired and permanently maintained procedure will prevent you from omitting some parts of the teeth. Do the cleaning for at least 2 minutes which is the minimal recommended time by stomathologists. You may use the timer **C1** for this purpose. If the plaque is removed properly, the gums will gain pink colour and will surround the teeth fast.

Timer (C1)

The timer in the form of hourglass allows to control the length of teeth cleaning (approx. 2 min.).

IV. MAINTENANCE

Clean the toothbrush regularly after each use! Do not use any rough and abrasive cleaning agents! Wipe the surface with a soft damp cloth. Rinse the brush head thoroughly for a few seconds under running water while leaving the propulsion unit on (pic. 3). Turn the propulsion unit by pressing the 0 button, take off the brush head and rinse both the parts in running water separately under running water (pic. 3) and let dry. To achieve hygiene cleanliness we recommend immersing the brush head every 2 - 4 in a cleaning agent for dentures.

Replacing the brush head

To achieve the best results while cleaning teeth, it is necessary to replace the brush head every time its bristles start to bend outward, i.e. after approx. 3 months. You can buy the replacement head at your distributor.

V. ENVIRONMENT

If dimensions permit, all pieces have printed symbols of materials used for production of packing, components and accessories as well as recycling symbols. The symbols on the product or documentation means the electric or electronic parts used in the appliance may not be disposed off together with communal waste. Hand over the electric or electronic part to a waste collection yard for free disposal. Proper disposal of the product will help preserving valuable natural resources and preventing potential negative impacts on environment and human health that would otherwise result from improper waste disposal. For more details, please contact your local authorities or collection yard. Fines may be imposed for improper disposal of this type of waste pursuant to national directives. Remove flat batteries from the appliance and dispose them off in a suitable way using special collection centres. Do not incinerate the batteries!

Specialized service is required for overhaul maintenance or maintenance requiring disassembly of the internal components of the appliance!

The right for warranty repair expires if the above mentioned instructions are not adhered to.

VI. TECHNICAL DATA

Weight (kg) about 0,1

Regarding EMC, the product conforms to the Council Directive No. 2004/108/EC as amended.

The manufacturer reserves its right for unimportant deviations from the standard design without influence of the product function.

DO NOT IMMERSE IN WATER.

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

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