SOEHNLE

Shape Sense Connect 200

D	BEDIENUNGSANLEITUNG	10
GB	USER MANUAL	14
NL	GEBRUIKSAANWIJZING	16
I	ISTRUZIONI PER L'USO	19
F	MODE D'EMPLOI	22
E	INSTRUCCIONES DE MANEJO	25
P	MANUAL DE INSTRUÇÕES	28
DK	BRUGSANVISNING	31
S	BRUKSANVISNING	34
NO	BRUKSANVISNING	37
FI	KÄYTTÖOHJE	40
CZ	NÁVOD K POUŽITÍ	43
SK	NÁVOD NA OBSLUHU	46
BG	УПЪТВАНЕ ЗА УПОТРЕБА	49
RUS	ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ	52
PL	INSTRUKCJA OBSŁUGI	55
ARAB	58 دليل الاستخدام	
EST	KASUTUSJUHEND	61
LV	LIETOŠANAS INSTRUKCIJA	64
LT	NAUDOJIMO INSTRUKCIJA	67
RO	INSTRUCȚIUNI DE UTILIZARE	70
SLO	NAVODILA ZA UPORABO	73
SR	UPUTSTVO ZA UPOTREBU	76
HR	UPUTE ZA UPORABU	79

D Interpretation der Messwerte • GB Interpretation of measured values • NL Interpretatie van de meetwaarden • I Interpretazione dei valori • F Interprétation des valeurs de mesure • E Interpretación de los valores medidos • P Interpretação dos valores medidos • DK Tolkning af målingerne • S Interpretation av mätvärden • NO Tolkning av måleverdiene • FI Mittausarvojen tulkinta • CZ Interpretácia výsledkov pomiaru • SK Interpretácia meraných hodnôt • BG Интерпретиране на измерените стойности • RU Интерпретация измерений • PL Interpretacja namierzonych wartości • ARAB شرح قيم القياسات • EST Väärtuste tähendus • LV Mērījumu skaidrojums • LT Matavimų vertčių paaiškinimas • RO Interpretarea valorilor măsurate • SLO Razlaga merilnih vrednosti • SR Interpretacija mernih vrednosti • HR Tumačenje mjernih vrijednosti

Display	1	2	3	4	5	6	7	8	9
	—			☺			+		
BMI	< 17,4	17,5 - 17,9	18,0 - 18,4	18,5 - 20,4	20,5 - 22,9	23,0 - 24,9	25,0 - 27,4	27,5 - 29,9	≥ 30
☼	< -12 %	- 12 %	- 7 %	-3 %	0 %	+ 3 %	+ 7 %	+ 12 %	> + 12 %
💧	< - 6 %	- 6 %	- 4 %	-2 %	0 %	+ 2 %	+ 4 %	+ 6 %	> + 6 %
👤	< - 6 %	- 6 %	- 4 %	-2 %	0 %	+ 2 %	+ 4 %	+ 6 %	> + 6 %

Alter/Age	♀						♂					
	☼		💧		👤		☼		💧		👤	
	—	☺	+	++	☺	☺	—	☺	+	++	☺	☺
17-18	< 15	≥ 15-25	≥ 25-33	> 33	> 52,5	> 28	< 8	≥ 8-18	≥ 18-24	> 24	> 58,5	> 38
18-30	< 20	≥ 20-29	≥ 29-36	> 36	> 49,5	> 28	< 8	≥ 8-18	≥ 18-24	> 24	> 57	> 38
30-40	< 22	≥ 22-31	≥ 31-38	> 38	> 47,5	> 28	< 11	≥ 11-20	≥ 20-26	> 26	> 56	> 38
40-50	< 24	≥ 24-33	≥ 33-40	> 40	> 46,5	> 28	< 13	≥ 13-22	≥ 22-28	> 28	> 55	> 38
50-60	< 26	≥ 26-35	≥ 35-42	> 42	> 46	> 28	< 15	≥ 15-24	≥ 24-30	> 30	> 54	> 38
60+	< 28	≥ 28-37	≥ 37-47	> 47	> 45	> 28	< 17	≥ 17-26	≥ 26-34	> 34	> 53	> 38

D Alle Angaben in % • GB All figures in % • NL Alle waarden in % • I Tutte le indicazioni sono in % • P Toutes les données en % • E Todos los datos en % • P Todas as informações em % • DK Alle tal i % • S Alla uppgifter i % • NO Alle tall i % • FI Kaikki tiedot prosentteina • CZ Všechny údaje v % • SK Všetky údaje sú v % • BG Всички данни в % • RU Все показатели в % • PL Wszystkie údaje jsou v % • ARAB نسبة البيانات بالمائة المئوية • EST Kõik väärtused %-na • LV Visi liekumi norādīti % • LT Visi skaičiai pateikiami procentais • RO Toate datele în % • SLO Vsi podatki v % • SR Svi podaci su u % • HR Svi podaci u %

BMI

D Body-Mass-Index • GB Body Mass Index • NL Body-Mass-Index • I Indice di massa corporea (BMI) • F Indice de masse corporelle • E Índice de masa corporal • P Índice de massa corporal • DK Body Mass Index • S Body-Mass-Index • NO Kroppsmasseindeks • FI Painoindeksi • CZ Wskaźnik masy ciała • SK Index telesnej hmotnosti • BG Body Mass Index • RU Body-Mass-Index • PL Index telesnej hmotności • ARAB مؤشر كتلة الجسم • EST Kehamassiindeks • LV Ķermeņa masas indekss • LT kūno masės indeksas (KMI) • RO Indexul de masă corporală • SLO Indeks telesne mase • SR Body Mass Index • HR Indeks tjelesne mase

D Körperfettanteil • GB body fat content • NL lichaamsvetaandeel • I percentuale di grasso • F taux de graisse • E porcentaje de la grasa corporal • P percentagem de gordura corporal • DK fedtprocent • S andelen kroppsfett • NO Fettprosent • FI Kehon rasvaprosentti • CZ podíl tělesného tuku • SK zawartość tłuszczu w ciele • BG Процент на телесната мазнина • RU доля внутреннего жира • PL test zsiřtartalmát • ARAB نسبة الدهون في الجسم • EST Rasvaprotsent • LV Ķermeņa tauku procents • LT kūno riebių procentinė dalis • RO Coeficient de grăsime corporală • SLO Delež telesne maščobe • SR Udeo telesne masnoće • HR Udio tjelesne masti


D Körperwasseranteil • GB body water content • NL lichaamswateraandeel • I percentuale acqua corporei • F taux d'eau dans le corps • E porcentaje del agua corporal • P percentagem de água corporal • DK vandprocent • S andelen kroppsvatten • NO Vannprosent • FI Kehon vesipitoisuus • CZ podíl vody v těle • SK zawartość wody w ciele • BG Процент на телесната вода • RU доля тканевой жидкости • PL test wżtartalmát • ARAB نسبة الماء في الجسم • EST Veeprotsent • LV Ķermeņa ūdens procents • LT kūno vandens procentinis • RO Coeficient de apă în corp • SLO Delež vode v telesu • SR Udeo telesne vode • HR Udio vode u tijelu


D Muskelanteil • GB muscle proportion • F pourcentage de muscles • I massa muscolare • NL Spierenaandeel • E porción de los músculos • P percentagem de músculo • DK muskler andelen • S muskelandel • NO Muskelprosent • FI Lihasprosentti • CZ podíl svalů • SK masy mięśniowej • BG Процент на мускулите • RU доля мышечной • PL izom arány • ARAB نسبة العضلات • EST Lihasmassi osakaal • LV Muskuļu procentis • LT raumenų procentinė dalis • RO Coeficient de masă musculară • SLO Delež mišične mase • SR Udeo mišića • HR Udio mišića


A

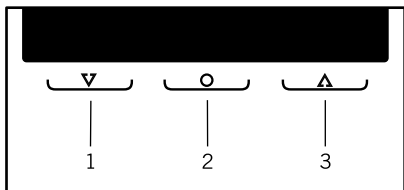
B

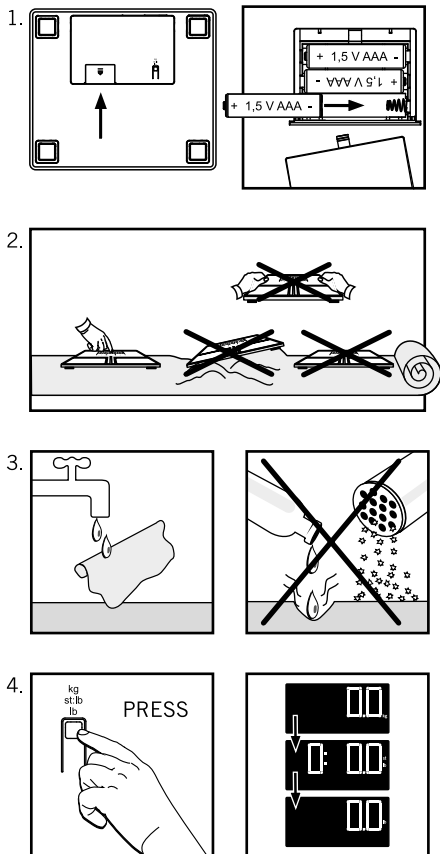
C

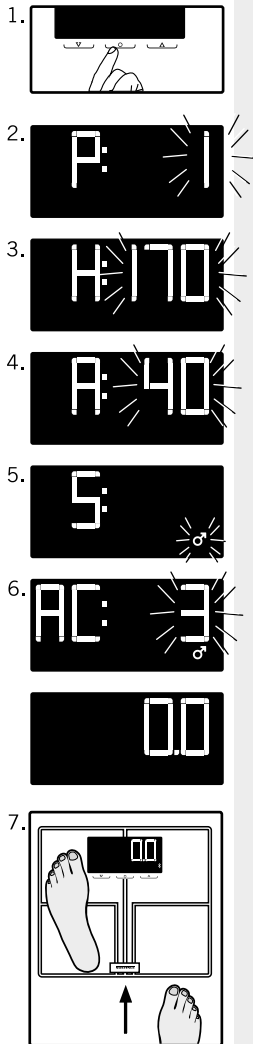D

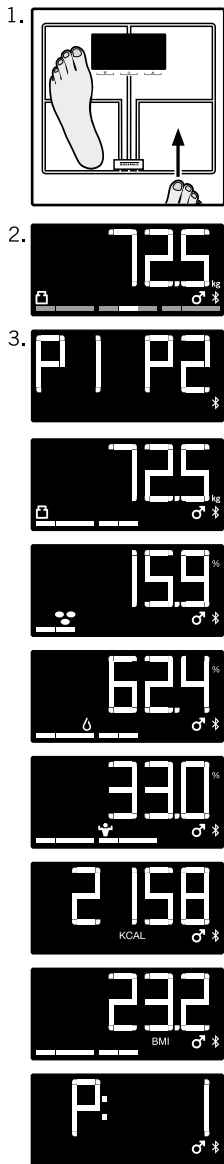E

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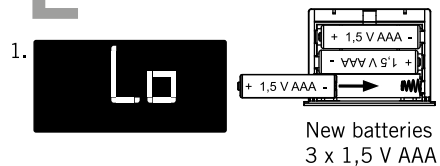GB This device may be used by children over the age of 8 and persons with reduced physical, sensory, or mental abilities or with a lack of experience and knowledge only if such are supervised or have been instructed in the safe use of the device as well as the resulting risks. Cleaning and maintenance shall not be made by children without supervision. Children shall not play with the device.

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Apple, iPhone, iPad, iOS, Apple Health and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Android, Google Fit, Google Play and the Google Play logo are trademarks of Google Inc.

The **Bluetooth®** word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

A. Operating elements

- 1. Minus (▼)
- 2. Confirm (○)
- 3. Plus (▲)

Sensor touch:
Only touch the keys softly!

B. Preparation

- 1. Insert batteries (3 x 1.5 V AAA).

Note: Immediately (as long as „0.0“ is still being displayed) place the body analysis scale on a level surface and wait until the scale powers down automatically.
Then begin with data entry. Otherwise, it might be possible that the weight of the scale would erroneously be added to the initial weighting process.

The same might occur during regular operation if the scale is held in the hand prior to the weighting process.

- 2. For all measurements, place scale level on a hard surface (not on carpet).
- 3. Cleaning and care: Clean only using a slightly damp cloth. Do not use solvents or abrasives. Never submerge scale in water.
- 4. A country-specific option to switch from kg/cm to st/in or lb/in is possible by pushing the key on the bottom of the scale to switch units

C. Data input

- 1. Press ○ button.
- 2. Use ▼/▲ buttons to select memory slot and confirm with ○ . → P
- 3. Use the ▼/▲ buttons to set your height and confirm with ○ . → H
- 4. In the next step, set your age and confirm with ○ . → A
- 5. Set gender and confirm with ○ . → S
- 6. To finish, set the activity level and confirm with ○ . (Hold down key for 3 sec.). → AC

AC	Explanation	Hrs per week
1	Hardly any physical activity	
2	Slight, active, sitting and standing activity	
3	House and garden work, mainly standing, sitting only occasionally	
4	Athletes and people with plenty of movement	Intensive training 5-10 hrs/week
5	Heavy manual labourer, highly active athletes	Intensive training at least 10 hrs/week

Athlete mode:

By setting the personal activity level, the respective calorie requirement can be precisely determined.

- 7. Now **immediately** place the body analysis scale on the floor, wait until the screen displays „0.0“ and then step on the scale with bare feet.

The initial measurement is indispensable and serves to store the data for later automatic recognition of persons. If no person steps on the scale, data input must be repeated. The first body analysis is performed after a person steps on the scale.

This scale is equipped with a Time-Out mode. It will automatically switch off after approx. 40 seconds if no keys are pressed.

Alternative: You can also perform this step conveniently via the Soehnle Connect app, in which case the data is transferred to the scale. The first measurement is necessary, however. Once you have entered the data on the scale for the recognition of persons as described under C, the Soehnle Connect app will automatically adopt this.

D. Body analysis

- 1. Place scale on a solid, level floor and wait until scale switches off automatically.
Step on the scale with bare feet. Stand still.
- 2. Body weight and the other calculated data are displayed.
- 3. If the determined values cannot be clearly allocated to one user, the possible memory slots will be displayed.

If the calculated weight can be allocated to more than 2 users, other possible memory slots will be displayed after 3 seconds.

Use the ▼/▲ keys (▼= left, ▲= right) to confirm the memory slot displayed there.

The scale measures the individual body values by means of bioimpedance analysis (body fat percentage, body water, muscle mass), calculates the basic metabolic rate per day (kcal) and the BMI (Body Mass Index). The 9-segment status bar shows the respective deviation from the normal value (see overleaf).

The scale switches off automatically.

Note: If needed, the user can step off the scale as soon as the first value appears on the display.

Interpretation of measured values

In order to be better able to assess your measured values, at the front of this operating manual you will find an overview table, from which you can take the commonly accepted optimal values.

The optimal range is marked in the table with a 😊. Values outside this range are either too low (→) or too high (+/+/+).

Here is an example:

A female user is 35 years old and is shown a body fat content of 28.2 % from the scale. For the age group of 30-40 year old women, this is a value in the optimal range (22-31 % → 😊).

Classification of weight according to BMI (source: World Health Organisation, 2016)

BMI (kg/m ²)	Classification
< 18.5	Underweight
18.5 - 24.9	Normal weight
25.0 - 29.9	Overweight
30.0 - 34.9	Obesity level I
35.0 - 39.9	Obesity level II
> 40.0	Obesity level III

E. Notifications

1. Replace battery (Lo = low battery)
2. Overloading of scale (max. 180 kg)
3. The measured value is above (H) or below (L) the threshold value specified below.

	L	🔴	H	L	🔴	H	L	🔴	H
♀	< 10 %		> 60 %	< 5 %		> 90 %	< 5 %		> 90 %
♂	< 4 %		> 60 %	< 5 %		> 90 %	< 5 %		> 90 %

F. Other add-on functions via the Soehnle Connect app

You can connect this body analysis scale with your smartphone (iOS® or Android™) via **Bluetooth®**. Use the Soehnle Connect app for this to make use of more functions among other things:

1. Difference to last weighing
2. Difference to the desired weight

G. Installation of the Soehnle Connect app

Install the Soehnle Connect app on your smartphone or tablet as follows:

Android™ smartphones / tablets

If you use an Android device, download and install the Soehnle Connect app from the Google Play Store. To do so, scan the QR code or open the app with the name "Play Store" and look for "Soehnle Connect app". Then load and install the free app.

iPhone® / iPad® (iOS®)

If you use an Apple iPhone or iPad, download and install the Soehnle Connect app from the Apple App Store. To do so, scan the QR code or open the app with the name "App Store" on your iPhone and look for "Soehnle Connect app". Then load and install the free app.

System requirements for using the Soehnle Connect app and full functionality:

Android™: 6.0 or later

iOS®: iOS 9 or later

Bluetooth®: 4.0 or later

Compatibility

The Soehnle Connect app is compatible with Google Fit® and Apple Health®.

You can find more information on the Internet at www.soehnle.de. Or contact our consumer service.

H. Connect scale with smartphone / tablet

1. Activate **Bluetooth®** on your smartphone or tablet and start the Soehnle Connect app.
- 2.a. When you start the Soehnle Connect app for the first time, you are automatically guided through the process for connecting a Soehnle appliance. Make sure that the scale is switched on (step onto the scale briefly).
- 2.b. Once you have already used the Soehnle Connect app, go to the connection in the Soehnle Connect app on the "main menu" and press "my devices" and then also "Add device".

I. Technical data

Bluetooth®: 4.1

Frequency band: 2.4 GHz

Maximum transmission power: 10 mW

Max. 200 kg / 100 g

Max. 440 lb / 0.2 lb

Max. 31 st / 0.2 lb

Age: 17-99 years

Batteries required: 3 x 1.5 V AAA

Guarantee

Leifheit AG offers a 3 year guarantee for this quality product as of the purchase date (or in case of orders, as of receipt of the goods). Guarantee claims must be asserted within the guarantee period immediately after occurrence of the defect. The guarantee includes product quality.

The following are excluded from the guarantee:

- (1) Wear caused by use or other natural wear;
- (2) Damage caused by improper use or handling (e.g. impact, shock, drop),
- (3) Damage caused by non-compliance with these operating instructions;
- (4) Battery, or rechargeable battery

In case of a guarantee claim, Leifheit, at its own discretion, either repairs the defective components or replaces the product. If a repair cannot be performed and an identical product is no longer available in the product range for the purpose of replacement, you will receive a replacement product of the closest equivalent value. Reimbursement of the purchase price is not possible in case of a guarantee claim. This guarantee furthermore does not grant any claims for damages.

To assert the guarantee, please contact the seller where you have acquired the product and present the defective product and the purchase receipt (copy). This guarantee is valid worldwide.

Your statutory rights, in particular warranty rights, continue to apply and are not limited by this guarantee.

EC conformity

Soehnle herewith declares that this device is in compliance with the basic requirements and the other pertinent regulations of the Directive 2014/53/EU. You can find the EU declaration of conformity at www.soehnle.com.

Battery disposal, EC Directive 2008/12/EC

Batteries are not a part of your regular household waste. You must return batteries to your municipality's public collection or wherever batteries of the respective type are being sold.

Batteries that contain hazardous substances feature this sign:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

Disposal of electric and electronic devices, EC Directive 2002/96/EC

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

Consumer service

The following contact partners are happy to help you in case of questions and comments:

Germany (08 00) 5 34 34 34
Monday to Friday 08:30 to 12:00

International soehnle.com

Legal information

Apple, iPhone, iPad, iOS, Apple Health and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Android, Google Fit, Google Play and the Google Play logo are trademarks of Google Inc.

The **Bluetooth®** word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

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Un utilizator are 35 de ani și primește un procent de grăsime corporală de 28,2% de la cântar. Acesta este un interval optim pentru grupa de vârstă 30-40 de ani la femei (22-31 % → ©).
The **Bluetooth*** word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.
Android, Google Fit, Google Play i logotip trgovine Google Play su žigovi tvrtke Google Inc.
Verbalni žig **Bluetooth*** i logotipovi registrirani su žigovi u vlasništvu tvrtke Bluetooth SIG, Inc.
LEIFHEIT AG