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PEP 257

NÁVOD NA POUŽITIE

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INSTRUCTION MANUAL

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Gallet

Osobní váha s tělesnou analýzou
Osobná váha s telesnou analýzou
Waga łazienkowa z analizą ciała
Bathroom body analysis scale
Személyi mérleg test analízissel

CE

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Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.


GENERAL SAFETY INSTRUCTIONS

- The appliance is designed exclusively for private use and for the envisaged purpose.
- This appliance is not fit for commercial use.
- Keep it away from sources of heat, direct sunlight, humidity (never dip it into any liquid) and sharp edges. Do not use the appliance with wet hands.
- If you should leave the workplace, always switch the machine off.
- Do not try to repair the appliance on your own. Always contact an authorized technician.
- Pay careful attention to the following “Special Safety Instructions for this unit“.

SPECIAL SAFETY INSTRUCTIONS FOR THIS UNIT

- Operate the appliance only in accordance with the specifications on the rating plate.
- Do not put the appliance on the wet places.
- Clean the scale with a slightly damp cloth. Do not immerse the scale in the water or use chemical / abrasive agents.
- Always use the scale on a hard, flat surface. DO NOT use on carpet
- DO NOT use the scale on wet and slippery floor. DO NOT step on the edge of the scale. The scale platform is slippery when wet. Keep it dry!.
- Always weigh yourself with the same scales places on the same surface. Uneven floors may affect the reading.
- Weigh yourself without clothing and footwear, before meals and always at the same time of day.
- Pacemaker users / wearers, or other medical implants with electronic components, are not recommended to use this scale.
- The body fat & hydration% are for reference only. Please always consult your doctor for advice of medication or diet in order to achieve ideal body fat level.
- To make sure accuracy of testing, please DO NOT use the scale in environment with strong electro-magnetic field.
- DO NOT use the scale on wet and slippery floor. DO NOT step on the edge of the scale. The scale platform is slippery when wet. Keep it dry!
- DO NOT use this scale on a place with intense shake. Keep stand still on scale when testing.
- Do not drop the device, or jump on the device when using. Avoid strong impact and sudden shake as it is a precision measuring device. The construction is delicate and sensitive to all unwanted external impact.
- This appliance is not intended for use by persons (including children) with reduced physical. Sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.

DESCRIPTION OF THE CONTROLS

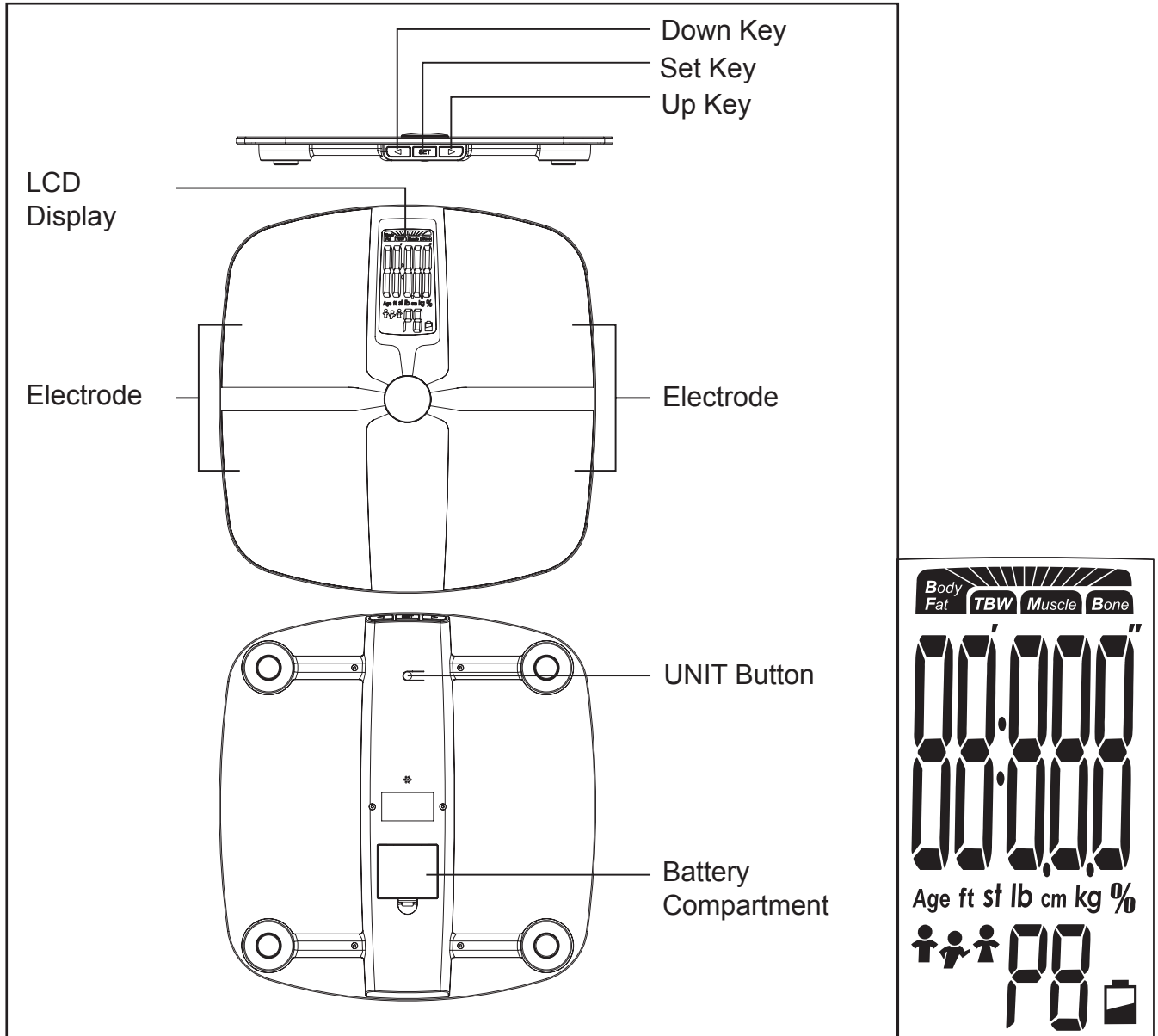
	Body Fat Analysis Result	kg	Kilogram
	Total Body Water Analysis Result	%	Percentage
	Muscle Mass Analysis Result		Man
	Bone Mass Analysis Result		Sportsman
Age	Age		Woman
ft	Foot		Sportswoman
st lb	Stone		User ID (Range from P1 to P4)
lb	Pound		Low Battery
cm	Centimeter		


Insert the Batteries

- Open the battery door in the back of the scale.
- Insert the batteries (4 x 1.5V AAA) into the battery compartment according to the polarity indications marked inside the compartment.

* The digits “88888” will be shown on the LCD.

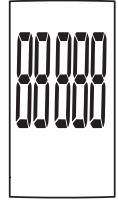- Close the battery door and wait until the digits “0.0kg” are shown on the LCD.

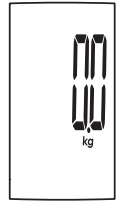CAUTION


- When the symbol  appears, the device will power off in four seconds. Then you shall replace with a new set of batteries. Please replace all four batteries in the same time. Do NOT mix the old batteries with the new one.
- Worn batteries are hazardous waste. Do NOT dispose of them together with the household garbage. Please refer to the local ordinances and recycling instructions regarding disposal of the worn batteries.

Set Up Your Profile

The glass body fat analyzer supports multiple users (Up to Four). You may follow below instructions to assign User ID and set up your own profile, including Gender, Stature, and Age.


1. Assigning User ID

- With batteries correctly installed, press “**SET**” key to enter setting.
- The system will request User ID selection first. As pictured below, “P1” blinks. The operator may press the function key ◀ or ▶ to select User ID among P1, P2, P3 and P4.
- Press “**SET**” key to confirm User ID.

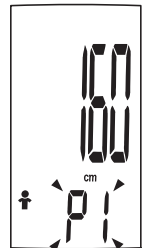
2. Setting Gender

- After confirming User ID, the system will divert to Gender setting automatically.
- As pictured below, the portrait  blinks. The operator may press the function key ◀ or ▶ to select Gender and Status.
- Press “**SET**” key to confirm User ID.

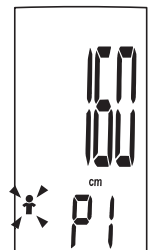
3. Setting Stature

- After confirming Gender, the system will divert to Stature setting automatically.
- As pictured below, the digits “160” blinks. The operator may press the function key ◀ or ▶ to increase or decrease the numeral. (Range: 100 cm to 220 cm).
- You may press and hold the function key ◀ or ▶ for fast changing the numeral.
- Press “**SET**” key to confirm stature.

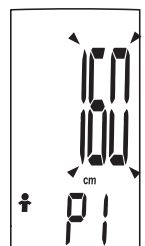4. Setting Age


- After confirming Stature, the system will divert to Age setting automatically.
- As pictured below, the digits “30” blinks. The operator may press the function key ◀ or ▶ to increase or decrease the numeral. (Range: 10 to 80 years old)
- You may press and hold the function key ◀ or ▶ for fast changing the numeral.
- Press “SET” key to confirm Age.

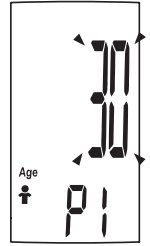When the LCD displays “0.0kg”, you may start measuring.

Select Measurement Unit


With batteries correctly installed, press “UNIT” button in the back of the scale to select measurement unit. The default measurement unit is “kilogram”. You may press “UNIT” button to choose among kilogram, stone and pound.

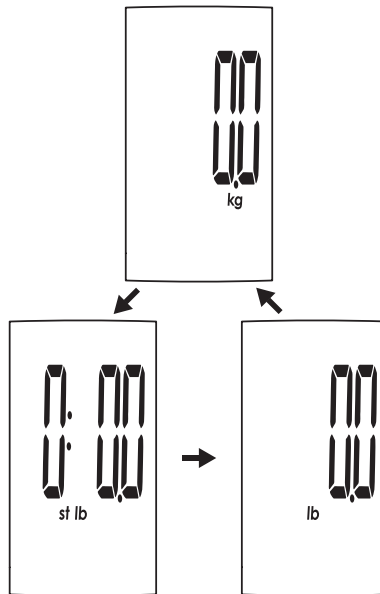
First Measurement

- Please process the user setting before your first measurement.
(Refer to Set Up Your Profile for more details)

STEP 1: Step on the platform barefooted.

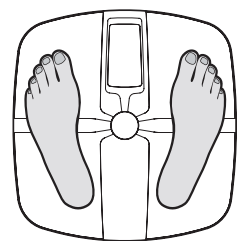
STEP 2: Stand still and keep full contact with the electrodes until the progress bar stops flashing and the measuring results displays.

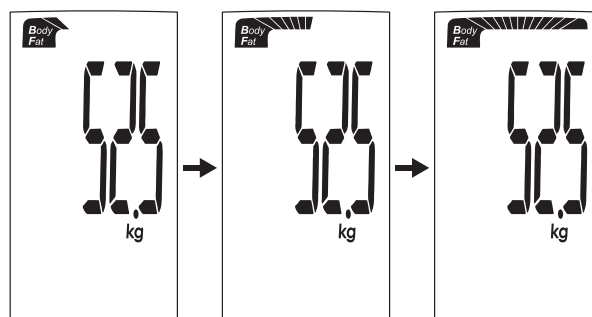STEP 3: Then the measuring results will be displayed sequentially for three times: Body Fat,


Total Body Water, Muscle Mass, Bone Mass.

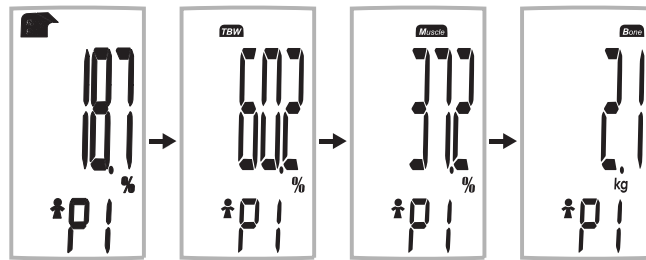
* If it fails to complete the analysis, the LCD will only display the weight data.

Daily Measurement


- With original SENSE ON patent technology, appliance will automatically switch on as you step on the platform barefooted.
- Stand still and keep full contact with the electrodes until the progress bar stops flashing.

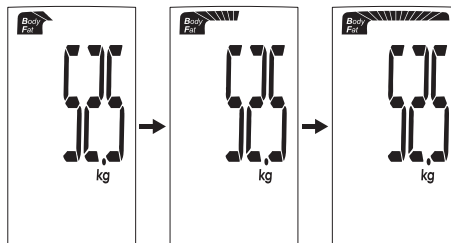- According to the analysis results, the system will automatically identify the possible User ID with most similar history records. Then the measuring results will be displayed sequentially for three times.

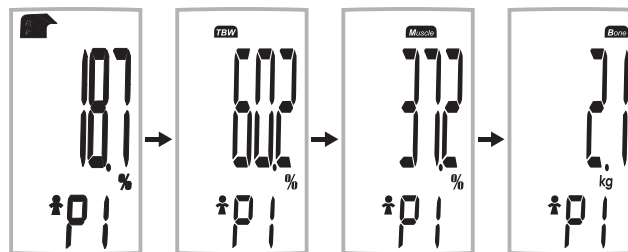
- When the system finds out two or more users with similar history records, it will notify you to choose between, for example, P1 and P2. You may press ◀ key for P1 and or ▶ key for P2. The measuring results will then be sorted into the User ID you selected and displayed sequentially for three times.

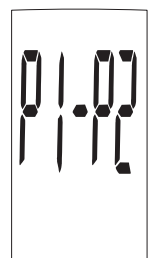
- If it fails to identify the possible User ID, the LCD will only display the weight data. (To find out the solutions, please refer to Troubleshooting for more details).

Body Fat Indicator

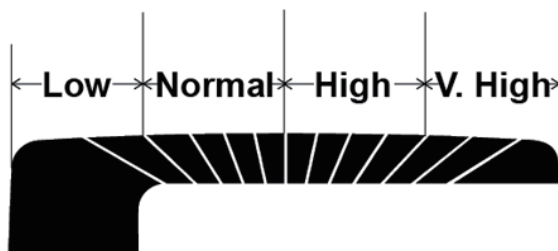The Progress Bar at the top of the LCD Display also functions as the indicator of body fat level. When the measuring result is displayed after analysis, the varying length of the bar indicates different levels of body fat.

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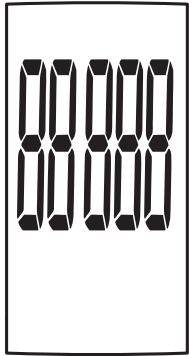Error	Description	Solution
	Overload. The device will power off in four seconds.	Stop using this scale for measurement.
	Low Battery. The device will power off in four seconds.	Replace all four batteries in the same time. Please purchase the authorized batteries for replacement.
	Install Batteries.	Install the batteries, the LCD can display and the device can start measuring.

When Measuring ...

Problem	Root Cause	Solution
Abnormal measuring results: - Too high; OR - Too low; OR - Huge difference between two recent measurement.	Incorrect posture	Please step on the platform barefooted and stand still.
	The device is located on the soft ground such as a carpet OR on a rugged surface.	Please place the device on a flat, hard surface.
	Cold body that may results in bad blood circulation.	Warm up your hands and feet to resume blood circulation and then measure again.
	Cold Electrodes.	Place the device in a warm room for a while and then measure again.
	Either your hands or your feet are too dry.	Wipe your feet with a damp cloth, keeping them slightly damp when starting measurement.
No display on LCD when the device powers on.	Batteries not yet installed.	Install the batteries. (Please refer to Insert the Batteries)
	Worn batteries.	Replace all four batteries in the same time. Please purchase the authorized batteries for replacement.
CANNOT proceed to analyze body fat, total body water, muscle mass and bone mass.	Step onto the platform wearing socks or shoes.	Please keep barefooted during the measurement, and keep full contact with the electrodes as well.
	The system cannot identify the possible User ID with most similar data.	Please assign a User ID following the instruction in Set Up Your Profile.
	The user fails to select the User ID from what the system found.	Please assign a User ID following the instruction in Set Up Your Profile.
The device powers off automatically.	Low battery.	Replace all four batteries in the same time. Please purchase the authorized batteries for replacement.

Table of Body Fat Level (Unit: %)

Age	Man / Sportsman				Woman / Sportswoman			
	Low	Normal	High	V. High	Low	Normal	High	V. High
20-29	<13	13,1 – 20	20,1 – 23	>23	<19	19,1 – 28	28,1 – 31	>31
30-39	<14	14,1 – 21	21,1 – 24	>24	<20	21,1 – 29	29,1 – 32	>32
40-49	<16	16,1 – 23	23,1 – 26	>26	<22	22,1 – 30	30,1 – 33	>33
50-59	<17	17,1 – 24	24,1 – 27	>27	<23	23,1-31	31,1 – 34	>34
60+	<18	18,1-25	25,1-28	>28	<24	24,1 – 32	32,1 – 35	>35

Table of Body Water Level (Unit: %)

Gender	Body Fat Percentage Range	Optimal Total Body Water Percentage Range
Men	4-14	70-63
	15-24	63-57
	22-24	57-55
	≥25	55-37
Women	4-20	70-58
	21-29	58-52
	30-32	52-49
	≥33	49-37

Bone Weight Chart

Weight	Female			Male		
	< 45 kg	45-60 kg	> 60 kg	< 60 kg	60-75 kg	> 75 kg
Ref. bone weight	~ 1,8 kg	~ 2,2 kg	~ 2,5 kg	~ 2,5 kg	~ 2,9 kg	~ 3,2 kg

The body muscle ratio(%): 25%-75%

Age	Female	Male
	Healthy	Healthy
10~99	>34	>40

CLEANING AND MAINTENANCE


- 1) Use a dry, soft fabric to clean the device. Or a lightly damp soft fabric with mild detergent to wipe clean the device. Do not use strong chemical solution such as alcohol, benzene for cleaning.
- 2) It could be slippery when the unit platform surface is wet, or semi-damp. Ensure the weighing platform surface is dry before using for safety reason.

TECHNICAL SPECIFICATIONS

Measurement Unit	Kilogram / Stone/ Pound
Measurement Range	5kg to 180kg / 0st: 11lb to 28: 5lb / 11lb to 397lb
Division	0.1 kg / 0.2 lb
Accuracy	0-50kg: $\pm 0.3\text{kg}$; 100-150kg: $\pm 0.5\text{kg}$; 50-100kg: $\pm 0.4\text{kg}$; 150-180kg: $\pm 0.7\text{kg}$
Working Environment	Temperature: 0 °C to 40 °C Humidity: $\leq 90\%$ RV
Storage Environment	Temperature: -20 °C to 60 °C Humidity: 10%RV to 93% RV
Power Source	6V (Four AAA-size Alkaline Batteries)
Auto-ON	SENSE ON technology
Auto-OFF	8 to 12 seconds after clearing Otherwise, 13-17 seconds after last operation

We reserve the right to change technical specifications.

Old electrical appliances, used batteries and accumulators disposal

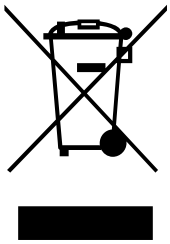This symbol appearing on the product, on the product accessories or on the product packing means that the product must not be disposed as household waste. When the product/ battery durability is over, please, deliver the product or battery (if it is enclosed) to the respective collection point, where the electrical appliances or batteries will be recycled. The places, where the used electrical appliances are collected, exist in the European Union and in other European countries as well. By proper disposal of the product you can prevent possible negative impact on environment and human health, which might otherwise occur as a consequence of improper manipulation with the product or battery/ accumulator. Recycling of materials contributes to protection of natural resources. Therefore, please, do not throw the old electrical appliances and batteries/ accumulators in the household waste. Information, where it is possible to leave the old electrical appliances for free, is provided at your local authority, at the store where you have bought the product. Information, where you can leave the batteries and accumulators for free, is provided to you at the store, at your local authority.

AUTO STEP ON – just step on the scale and stand still, the scale will automatically turn ON and start taking measurements.

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